


Accessories

Locking Dialysis Membrane Clamps

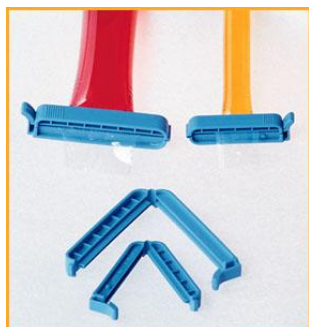

Superior in sealing dialysis tubing, our clamps stay locked until you release them, and they do not float. Constructed of Polyamide 6.6 from DuPont (NYLON). Available in three sizes to accommodate flat widths up to 46 mm, 65 mm or 105 mm. Ten clamps per package.

Part No.	Max Sealing Width	Total Length
CB-10100	105 mm	130 mm
CB-1050	46 mm	65 mm
CB-1070	65 mm	90 mm

Filter Funnels



Make loading and unloading dialysis tubing easy: Cut tubular membrane to length and clamp at one end, slip the open end onto the polypropylene funnel and secure with the silicone loop fastener. Open snap cap, fill with sample, and hang in container. Funnels are supported by a three-notched universal hanger that fits over the edge of common vessels.

Part No.	Part	Description
MFPS-05	Mini Filter Funnels	For 10 mm and 19 mm flat widths. Package of five funnels and two silicone loops.
MFPL-05	Macro Filter Funnels	For 23 mm, 25 mm, 32 mm, 40 mm, and 45 mm flat widths. Package of five funnels and two silicone loops.



Dialysis Tubing Bag Holders

These strong polypropylene holders suspend dialysis bags in containers. Easy to mark on top for identifying sequence when dialyzing multiple samples. Shown above with Locking Dialysis Membrane Clamp (not included). Package of 5 holders.

Part No.

MBHS-05

Accessory Kit with Membranes



Our dialysis start-up kit comes with a sampling of our handy accessories and two types of Cellu-Sep Regenerated Cellulose Tubular Membrane.

Each kit includes:

- 2 meters of Cellu-Sep T2, Nominal MWCO 6,000-8,000, 10 mm flat width.
- 2 meters of Cellu-Sep T3, Nominal MWCO 12,000-14,000, 25 mm flat width.
- Six Locking Dialysis Membrane Clamps, 46 mm max sealing width, 65 mm total length.
- 1 Dialysis Tubing Bag Holder.
- 1 Macro Filler Funnel with silicone loop fastener. For use with 25 mm flat width dialysis tubing.
- 1 Mini Filler Funnel with silicone loop fastener. For use with 10 mm flat width dialysis tubing.

Part No.

DMAK-01

Cellu-Sep Flat Sheet Membrane Squares for the QuixSep 5ml Dialyzer and MicroDialyzers Available in three popular Molecular Weight Cut-Offs (MWCO):

- F1 – 3,500 MWCO
- F2 – 6,000 – 8,000 MWCO
- F3 – 12,000 – 14,000 MWCO

The sheets are packaged 10 to a box, and measure 55mm x 55mm. This provides 10 dialysis procedures for the QuixSep 5ml dialyzer... OR... if you are using the .1 , .5 and 1ml MicroDialyzers, each sheet can be cut into 4 equal pieces, providing up to 40 procedures per box.

Preparation is easy... just a brief rinse in distilled water is all that is necessary for most applications.

Part No.**Description**

Q5010-55	F1 Membrane, 3.5K MWCO 55mm x 55mm, 10/pk
Q8010-55	F2 Membrane, 6-8K MWCO, 55mm x 55mm, 10/pk
Q1210-55	F3 Membrane, 12-14K MWCO, 55mm x 55mm, 10/pk