


# accessories

## gel rulers

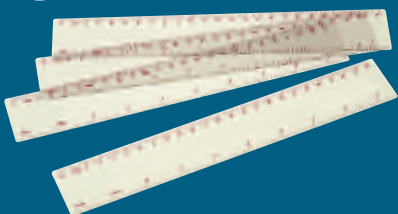

- SHARP, CLEAR, EASILY VISIBLE MARKINGS
- FOR USE WITH ETHIDIUM BROMIDE GELS
- IDEAL FOR GEL PRESENTATION


- COLOUR CODED FOR THICKNESS
- COMB AND GUIDE IN ONE PIECE FOR CONVENIENCE
- ADJUSTABLE HEIGHT

- EFFECTIVE COOLING
- COOLS WITHOUT EXTERNAL CHILLING SYSTEMS
- IMPROVES GEL RESOLUTION AND REDUCES RUN TIMES

## HORIZONTAL GEL SYSTEM


### buffer saver blocks

**Buffer saver blocks** physically reduce the volume in positive and negative compartments of a gel chamber and so reduces buffer requirements, saving cost. Blocks are available for all multiSUB Horizontal Gel Units.

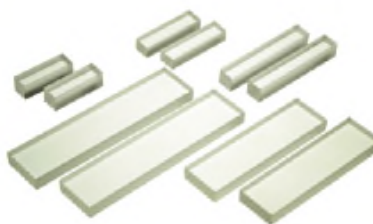


### gel levelling table


This **Gel Levelling Table** provides a useful way of pouring gels of uniform thickness. The table features a large 32 x 26cm platform with large thumb-wheel levelling feet at each corner. A levelling bubble is supplied.

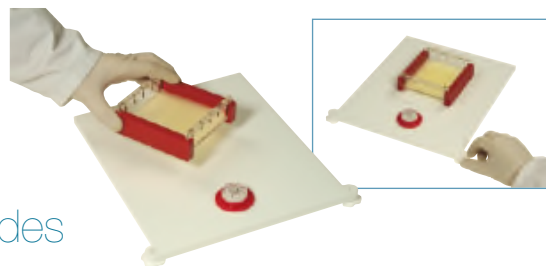


### combi-comb/loading guides

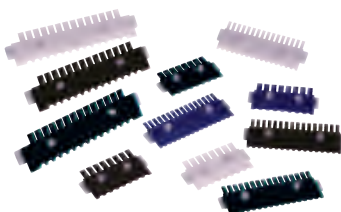

These **guides** combine the convenience of a regular comb with the appropriate guide (same teeth number and width) for loading samples, simply flip over.

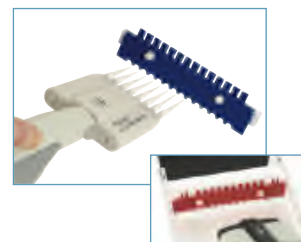

### high efficiency horizontal gel cooling

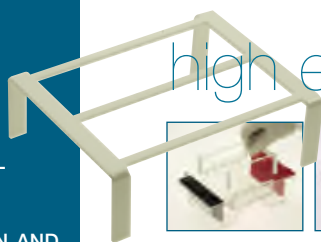

**Comprising a removable Gel Platform and CoolPack** this unique system will cool the buffer chamber of multiSUB horizontal gel units without any external chilling. Simply place the platform over the running tray and place a pre-chilled CoolPack on top. This allows the CoolPack to sit in the buffer, above the gel directly, removing any heat generated.

#### ORDERING INFORMATION

##### High Efficiency Horizontal Gel Cooling Systems

MSMINICP	CoolPack and Platform for MSMINI
MSMIDICP	CoolPack and Platform for MSMIDI
MSCHOICECP	CoolPack and Platform for MSCHOICE
MSMAXICP	CoolPack and Platform for MSMAXI
MSSCREENCP	CoolPack and Platform for MSSCREEN

##### Horizontal Buffer Saver Blocks

MSMINIBSB	Buffer Saver Blocks for MSMINI, pk/2 saves 100ml of buffer
MSMIDIBSB	Buffer Saver Blocks for MSMIDI, pk/2 saves 100ml of buffer
MSCHOICEBSB	Buffer Saver Blocks for MSCHOICE pk/2 saves 190ml of buffer
MSMAXIBSB	Buffer Saver Blocks for MSMAXI pk/2 saves 450ml of buffer
MSSCREENBSB	Buffer Saver Blocks for MSSCREEN, pk/2 saves 625ml of buffer

##### Horizontal Systems Combi Comb / Loading Guides

(x = comb well number, y = comb thickness)	
MS7-x -yLG	Combination Comb and Loading Guide for MSMINI
MS10-x -yLG	Combination Comb and Loading Guide for MSMIDI
MS15-x -yLG	Combination Comb and Loading Guide for MSCHOICE
MS20-x -yLG	Combination Comb and Loading Guide for MSMAXI
MS26-x -yLG	Combination Comb and Loading Guide for MSSCREEN

##### Gel Levelling Table

GLT	Gel Levelling Table
-----	---------------------

##### UV Fluorescent Ruler

RLR	UV Fluorescent Ruler
-----	----------------------