



- Repeatable Dial Speed Control
- Simple Operation
- Near Silent Whilst Running
- 5-50 RPM Speed Adjustable
- 350 x 350mm Platform
- Large Range of Accessories Available

Excellent Value with Unrivalled Reliability

The ERPM4 is the perfect everyday rocking platform for a wide range of mixing applications. The solid state speed control has a speed range of 5 - 50 RPM and is controlled from a simple to operate turn dial on the front of the mixer, whilst the Ratek silent belt drive system will deliver years of quiet reliable operation without service. The platform is large enough to accommodate up to 16 x 250mL flasks or 5kgs, yet small enough to be part of your daily workstation.



SPECIFICATIONS

ERPM4	
Mixing action	Rocking 7° to 15° (factory preset at 10°)
Platform size	350 x 350 mm
Speed range	5-50 RPM
Speed control type	Solid state with optical feedback
Motor	Long life super-quiet brushless AC
Transmission type	Silent Belt
Safety features	<ul style="list-style-type: none"> • Motion interrupt stop • Slow start • Unbalanced load stop
Power input	240V AC/40W
Overall dimensions	W370 x D370 x H200 mm
Nett weight	10.1 kg

Slow acceleration is incorporated into the drive system to reduce the risk of spillage, as well as a motion interrupt safety cut-out. The supplied non slip mat is ideal for shaking flat-bottomed containers, and a wide range of accessories is available for higher speed mixing or to retain more specialised containers. These include Ratek's amazing AT5 adhesive mat tray, TCT5 tulip clip tray and the RR5 universal rack.

The ERPM4 is simple to operate yet offers precise repeatable speed control in a package that is built to last.

Backed by Ratek's unparalleled service, support and warranty, the ERPM4 is an excellent value proposition for any lab.


RELATED PRODUCTS

Product Code	Description
AT5	Adhesive mat tray to suit medium mixers
RR5	Universal rack to suit medium mixers
TCT5	Tulip clip tray to suit medium mixers (clips sold separately)
DTT5	Dual-tier tray to suit medium mixers
MR5	Mojonnier rack tray to suit medium mixers

YOU MAY ALSO LIKE

RPM5 - Digital Rocker

The RPM5 is a medium sized rocking platform with precision repeatable digital control plus timer, tacho and user presets.

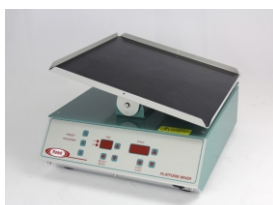

AUSTRALIAN
MADE & SUPPORTED

3 YEAR
PARTS WARRANTY



Australian distributors:
Fisher Biotec Australia
free call: 1800 066 077
email: info@fisherbiotec.com
web: www.fisherbiotec.com

